

Tutorial 5: Generating *in-silico* microscopy image with different hues

Subhamoy Mahajan

28 Dec, 2020

1. Generate *in-silico* monochrome image data files


The monochrome image data file generated in Tutorial 1 is used (“img100_lam670_fs800.dat” and “img100_lam518_fs800.dat”).

2. Generate colored *in-silico* microscopy images

First, the microscopy image is generated with red (0°) and green (120°) hues,


```
term$ python ../../mono2color.py -f img100 -p png_param_rg.dat -t -1
term$ mv img1000.png img_rg.png
```

For simplicity mono2color.py was run with “-t -1” option, so that the output is “img1000.png”, which is then renamed to “img_rg.png”.

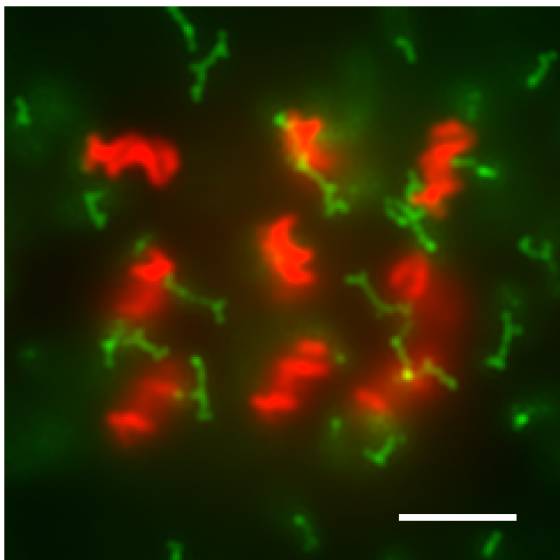

Second, the microscopy image is generated with orange (60°) and violet (270°) hues,

```
term$ python ../../mono2color.py -f img100 -p png_param_ov.dat -t -1
term$ mv img1000.png img_ov.png
```

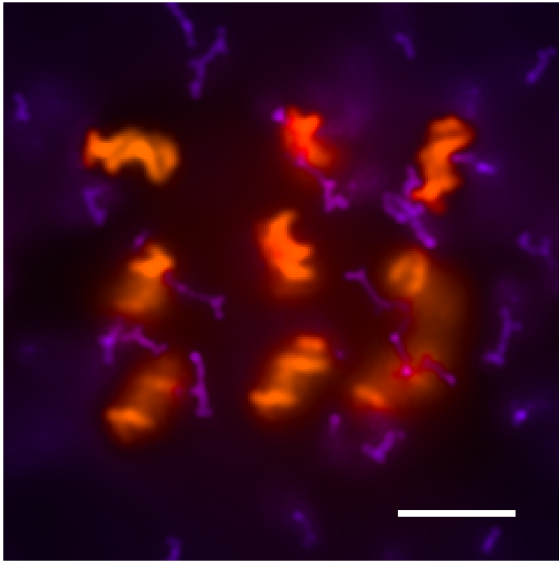

Third, the microscopy image is generated with cyan (180°) and magenta (300°) hues,

```
term$ python ../../mono2color.py -f img100 -p png_param_cm.dat -t -1  
term$ mv img1000.png img_cm.png
```

